belcils MED

Gotas oftálmicas INTENSIVE

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belcils MED

Eye drops INTENSIVE

Instructions for use - please read attentively

Composition

Sodium hyaluronate 0.15%, dexpanthenol 2%, sodium chloride, sodium citrate dihydrate, citric acid monohydrate, purified water.

Preservative-free. Phosphate-free.

Read the instructions for use.	MD Medical device.
Manufacturer.	Net content.
Expiry date.	Dispose of product at the Sigre point.
STERILE A Sterile, aseptic processing technique.	UDI Unique device identifier.
Caution, check the Instructions for use (for important warnings/precautions)	
LOT Batch number.	
CE mark.	
6 MO Expiry date after opening is 6 months	

Revision date: May 2025.

MD CE 1304

EXPIRY DATE: Printed on the packaging. Expiry date after opening is 6 months.
Do not use the product after its expiry date.

STORAGE: This medical device does not require any special storage condition.
Keep out of sight and reach of children.

PACKAGING: 10 ml

Distributed by:

Laboratorios Viñas
Provença 386 - 08025 Barcelona
ESPAÑA/SPAIN

JADRAN-GALENSKI LABORATORIJ d.d.
Svilno 20, 51000 Rijeka, Croatia

Belcils MED Eye drops INTENSIVE

PURPOSE

Belcils MED INTENSIVE eye drops relieve symptoms of mild to moderate dry eye. They stabilize the tear film, lubricate and soothe the eye surface.

The application of the eye drops provides intensive care and protection from:

· **Dryness and mechanical irritation due to external conditions**: wind, dry air, air-conditioning (AC), air pollution, dust, wearing lenses, diagnostic or surgical interventions on the eyes.

· **Eye strain and irritation due to lack of blinking** in long periods of focused viewing at: digital screens, driving, precise work.

Belcils MED INTENSIVE is double effective: it contains **sodium hyaluronate**, which effectively moisturizes the eye and stabilizes the tear film, and **dexpanthenol (provitamin B5)**, which additionally binds water, soothes and protects the eye surface.

The application of the drops provides fast relief of the unpleasant eye symptoms.

Dry eye symptoms can be relieved and thus complications can be prevented.

Risk factors for dry eyes, i.e., disorder of the tear film caused by lack of tears or their excessive evaporation are:

· Aging

· Gender (hormonal changes in women during pregnancy, breast-feeding, menopause)

· Modern lifestyle (working for long periods on a computer, living, working and travelling in air-conditioned spaces)

· Environmental conditions (wind, dry air, dust, smoke, air pollution, AC, neon lights)

· Eye diagnostic and surgical interventions

· Presence of a foreign object in the eye

· Wearing contact lenses

· Use of decongestive eye drops

· Medications (for treating depression, hypertension, for prevention of pregnancy)

· Diseases (rheumatoid arthritis, diabetes, asthma, thyroid problems, glaucoma)

WARNINGS AND PRECAUTIONS

Do not use Belcils MED INTENSIVE eye drops:

· In the case of hypersensitivity (allergy) to sodium hyaluronate, dexpanthenol or any other ingredient of Belcils MED INTENSIVE eye drops.

· In the case of an allergic reaction, stop using the product immediately and contact a doctor.

Other remarks:

· If you notice that the Belcils MED INTENSIVE bottle or applicator is damaged, do not use the product.

· To ensure safe application and hygiene, each bottle should be used by one person only.

· There are no known interactions between Belcils MED INTENSIVE and contact lens solutions, other eye drops or topical eye ointments. If you use another eye product, wait at least 15 minutes between applications.

Note:

Immediately after application of the drops, your sight might be blurred for very short period of time. It is therefore recommended to refrain from driving a car, operating a vehicle or machine and wait until your sight becomes sharp again.

PRODUCT SAFETY

Belcils MED INTENSIVE eye drops are preservative and phosphate free, which reduces the possibility of irritation and damage to the surface of the eye.

The dropper technology includes a filter which makes the product safe for use because it prevents the contamination of the eye drops and preserves sterility of the contents even 6 months after opening.

The product is generally well tolerated. However, if you notice any side effects, contact your doctor or pharmacist.

HOW TO USE BELCILS MED INTENSIVE

Always use Belcils MED INTENSIVE according to these instructions for use. Consult your doctor or pharmacist before use if you are not sure.

If not otherwise prescribed by your doctor, instil 1-2 drops in the eyes, several times a day, as needed.

Blink several times to distribute the drop evenly over the eye surface.

Suitable for long-term usage

Avoid contact of the dropper tip with the surface of the eye or any other surface.


Before the first application, it is necessary to activate the bottle. For this purpose, remove the stopper from the pump, turn the bottle upside down and **press the applicator several times until the first drop appears.**

Information for contact lens users


The product can be used while wearing contact lenses.

Instructions for use


1. Wash your hands thoroughly.
2. Remove the stopper from the pump (image 1).
3. Grasp the pump by hand as shown in image 2.
4. Turn the bottle upside down and press the applicator until the first drop appears (image 3).
5. With your head back, carefully instill 1 drop in the eye (image 4).
6. Replace the stopper on the pump.

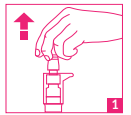1. Remove the stopper from the pump.

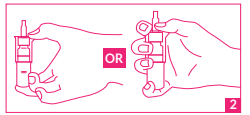2. Grasp the pump by hand in one of the two ways shown in Image 2

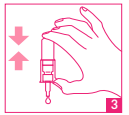3. Turn the bottle upside down and activate the pump by pressing it, until the first drop appears

4. Hold the bottle upside down and by pressing the pump, instil prescribed number of drops into the eye.

ASUDXXX-V1 (XX)